
Apparaat en bediening

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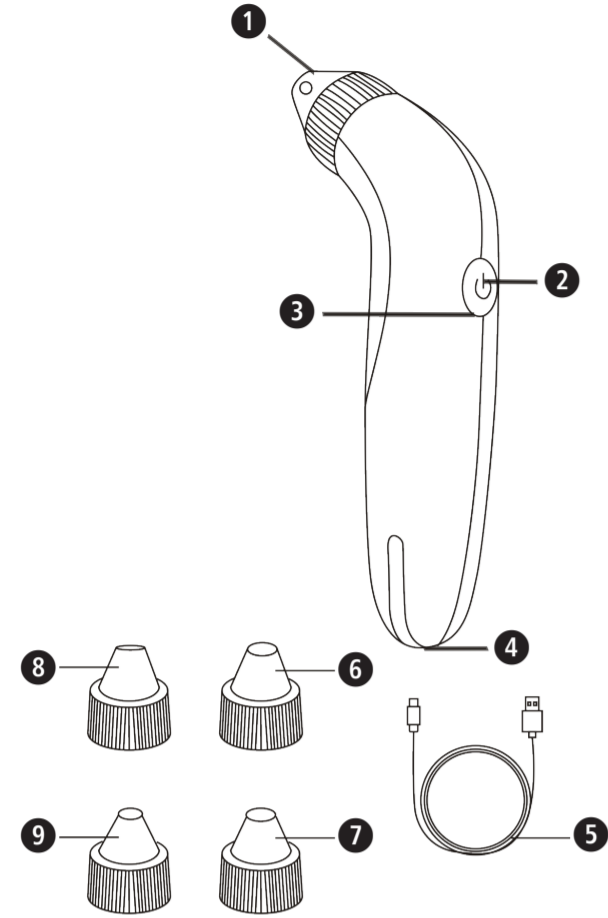
Jagenbergstr. 19
41468 NEUSS
DEUTSCHLAND

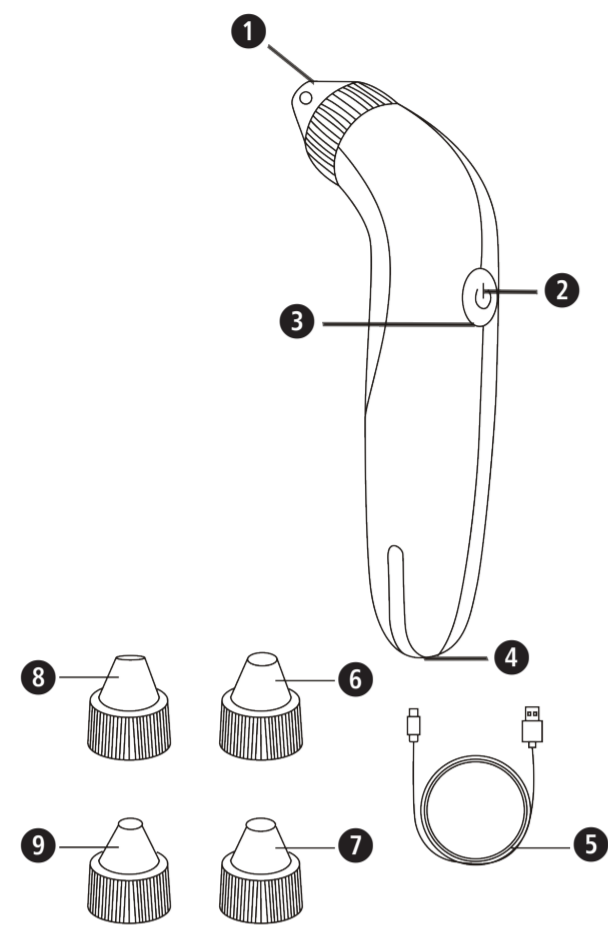Device and controls

- ➊ Position for exchangeable head
- ➋ On-/Off-button
- ➌ Charging indicator
- ➍ USB charging port
- ➎ USB charging cable
- ➏ Head (big, round opening)
- ➐ Head (thick and round)
- ➑ Head (elliptic, transparent)
- ➒ Head (small, round opening)

Explanation of symbols

IMPORTANT

Non-observance of these instructions can result in serious injury or damage to the device.

WARNING

These warning notes must be observed to prevent any injury to the user.

CAUTION

These notes must be observed to prevent any damage to the device.

NOTE

These notes give you useful additional information on the installation or operation.

Do not use the device in damp rooms (e.g. when bathing or showering).

LOT LOT number

Manufacturer

GB Safety instructions

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

- The device must only be used for its intended purpose as described in the instruction manual. Using the appliance for any other purpose invalidates the warranty.
- Do not immerse the device in water or other fluids and only use the device with dry hands.
- Always carefully inspect the device for damage before use. Do not attempt to use a defective unit.
- Do not dismantle or modify the device in any way.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Do not use this device near flammable or corrosive materials or liquids.
- Do not use this unit to supplement or replace medical treatment.
- Do not store this device in direct sunlight and store it only in dry and well-ventilated areas.
- If the supply cord of this appliance is damaged, it must only be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- The rechargeable battery must not be disassembled, heated, incinerated or exposed to water.
- During the use of this device no water or cream residues must come in contact with the device as this may lead to defects.
- Before use, remove any jewellery (e.g. necklaces, earrings etc.) and contact lenses.
- This device must not be used on sensible parts of the skin (e.g. eyes, ears, neck etc.) - risk of injuries!
- Do not use this device on parts of the skin, which have been altered by cosmetic procedures (e.g. after cosmetic surgery).
- Do not use the device during pregnancy or existing heart diseases.
- If you should feel unwell during use, stop the use immediately.
- Cleaning and head exchanging must only be done when the device is switched off.
- If the device has been fallen or other damages are present, do not use the device - contact the customer service!
- Do not charge the battery during the use!
- The device must not be used by patients with allergies, skin diseases, dermatitis, swellings, eczema or other serious (skin-) injuries.
- Cleaning and user maintenance shall not be made by children without supervision.

Scope of supply

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use the it and contact your dealer or your service centre.

The following parts are included:

- 1 **medisana** Pore cleanser **PR A60**
- 4 different heads
- 1 USB charging cable
- 1 Instruction manual

If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING

Please ensure that polythene packing is kept away from the reach of children! Risk of suffocation!

Charging the battery

Charge the battery for at least 1 hour before first use.

Explanation of the charging indicator ➌:

Red, constant light: Charging in progress - please do not interrupt.

Green light: Battery is fully charged.

Charge the battery in accordance with the following guidelines

- Charge the battery only with the supplied charging cable.
- Do not use the device during charging!
- **Charging the battery:**
Connect the mini USB plug of the included charging cable to the charging port ➍ of the device and connect the other end to a suitable power source.

Use

By pressing the On-/Off-button ➋ once, the device will be switched on at the lowest level (level 1). Pressing the button twice activates the medium level (level 2) and three presses activate the highest level (level 3). A fourth press on the button switches the device off again.

The head ➏ (big, round opening)

This head provides a strong suction for heavily contaminated pores. Keep moving slowly on the skin from bottom to top with activated device.

The head ➐ (thick and round)

This head is mainly suitable for the exfoliation of the facial skin. You should not use this head more than twice a week.

The head ➑ (elliptic, transparent)

This head is especially suitable for eignet sich insbesondere for wrinkle smoothing, promoting skin elasticity and blood circulation.

The head ➒ (small, round opening)

This head is designed for sensible skin. It features a less strong suction than the head ➏, but is also suitable for pore cleaning.

General instructions for handling:

- The application should not be longer than 30 seconds maximum twice daily at the same part of the skin.
- Hold the head with light pressure on the skin area only - do not apply force or greater pressure!
- After application, the skin area may become a little red. This is completely normal and disappears by itself after some time.
- After using the device, do not scrub the skin area with brushes or rags. You can, however, moisten the skin with a little moisturizer.

Cleaning and maintenance

- Before cleaning the unit, make sure it is switched off. Clean the unit with a soft, dry cloth.
- The head may be cleaned separately e.g. with a cotton swab.
- Under no circumstances should you use aggressive detergents, strong brushes, solvents or alcohol, as they may damage the surface of the unit.
- When cleaning, never immerse the unit in water and make sure that no water gets into the unit.
- You should ideally store the device in its original packaging and keep it in a clean, cool and dry place.

Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

Technical specifications

| | |
|--------------------|---|
| Name and model: | medisana Pore cleanser PR A60 |
| Battery: | Lithium-Ion, 3.7V, 400mAh (charging time approx. 1 h) |
| Input: | 5V ===, 0.4A |
| Power consumption: | 4 watts |
| Application time: | 45 minutes |
| Vacuum suction: | 60-65 kPa, level 1: 63.8 kPa ± 10%, level 2 : 64.1 kPa ± 10% kPa, level 3: 64.5 kPa ± 10% |
| Noise: | 58.5 dB (in 1 m distance) |
| Size approx.: | 19.5 x 8.5 x 4.5 cm |
| Weight approx.: | 170 g |
| Item number: | 99458 |
| EAN number: | 40 15588 99458 6 |

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The current version of this instruction manual can be found under www.medisana.com

Warranty and repair terms

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **medisana** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana GmbH

Jagenbergstr. 19
41468 NEUSS
GERMANY